
MARQUE: SOEHNLE

REFERENCE: 63806 EASY CONTROL

CODIC: 3575250



SOEHNLE

BODY BALANCE EASY CONTROL

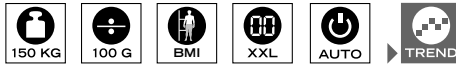BEDIENUNGSANLEITUNG
OPERATING INSTRUCTIONS
MODE D'EMPLOI
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING
INSTRUCCIONES DE MANEJO
MANUAL DE INSTRUÇÕES
BRUKSANVISNING
BRUGSANVISNING
KÄYTTÖOHJEET
KEZELÉSI ÚTMUTATÓ
INSTRUKCJA OBSŁUGI
NÁVOD K POUŽITÍ
ИНСТРУКЦИИ ПО ЭКСПЛУАТАЦИИ
KULLANMA KILAVUZU
Οδηγίες λειτουργίας
NAVODILA ZA UPORABO
UPUTA ZA UPORABU

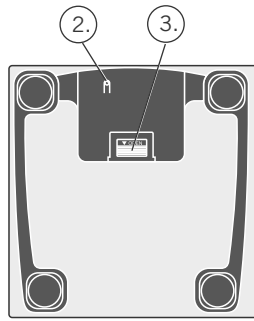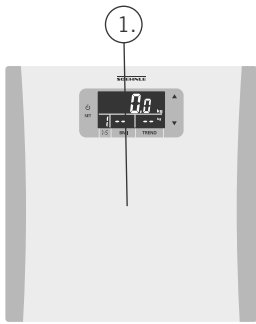DE

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EN

Start-up operation

1. Insert batteries (3 x 1.5 V AAA).
Note: stored data will be reset when batteries are replaced.
2. Set units.
3. Place scale on a hard, level surface.

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

EN

Entering personal data and weighing

1. Push ON/SET key.
Use ▲/▼ keys to select memory slot and confirm (SET).
2. Then use ▲/▼ keys to set body height and confirm (SET).
3. Finally enter the target

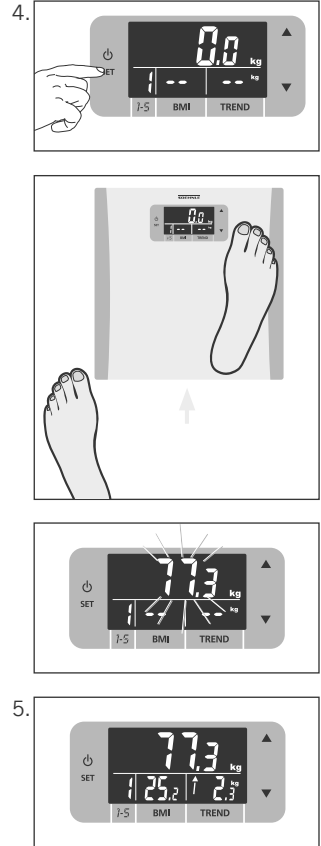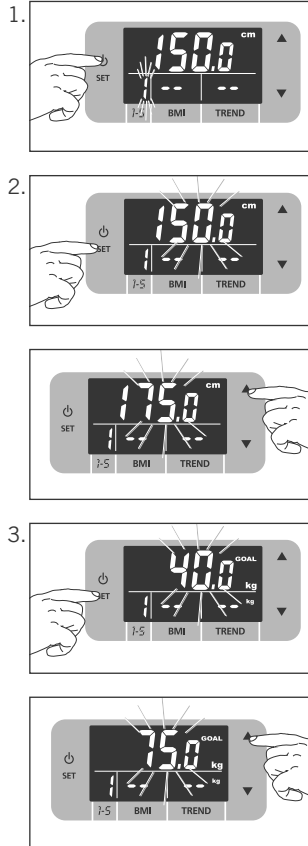- weight using the ▲/▼ keys and confirm (SET).
4. The 0-display will appear. Step on the scale for initial weighing. Please stand still for as long as the display flashes.
5. After determination of weight, the scale will display the body weight, the memory slot, the calculated body mass index and the

difference of your current to your target weight.

To change or delete data simply write over the existing data.

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EN

Weighing with BMI calculation and display of weight difference

1. Step on the scale.
2. The body weight will be determined and the automa-

tic recognition of persons will be performed. Please stand still for as long as the display flashes.

3. The scale will then display the body weight, the memory slot, the calculated body mass index and the difference of your current to your target weight.
4. If the body weight of other persons that are stored on the scale is approximately the same, a manual allocation will be required. This will be necessary if the memory slot number flashes. In that case simply select the memory slot using ▲/▼ keys and confirm (SET).

The scale will then display the body weight, the memory slot, the calculated body mass index and the difference of the current to the target weight of the correct person.

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Leifheit AG